
GROWNSY Electric Baby Nasal Aspirator

Model BC023

User Manual – Multilanguage Edition

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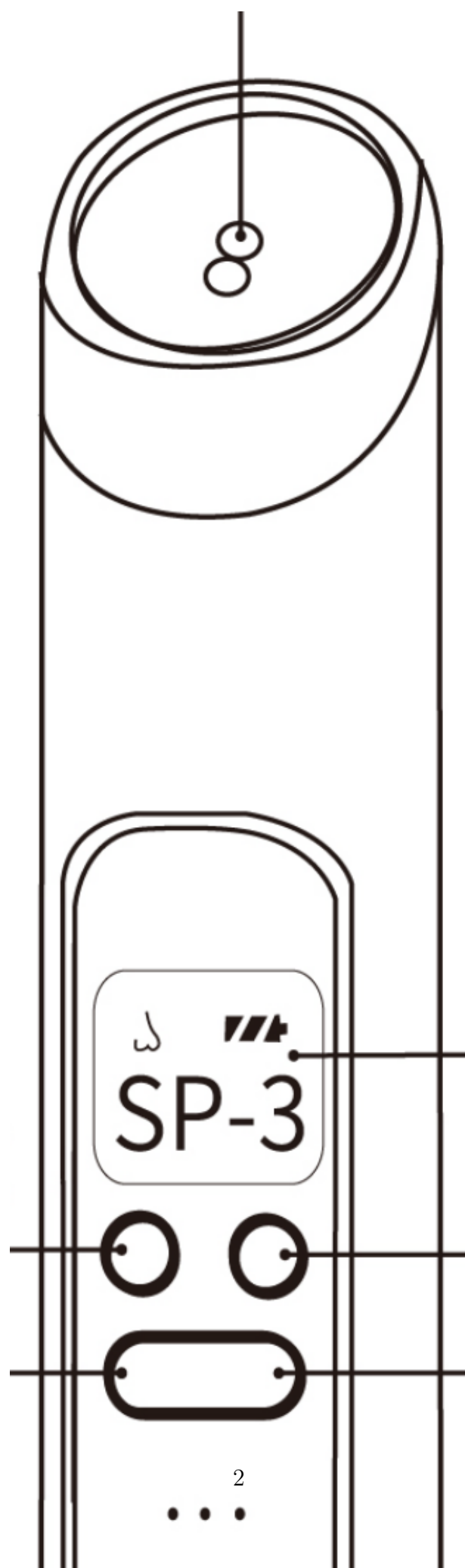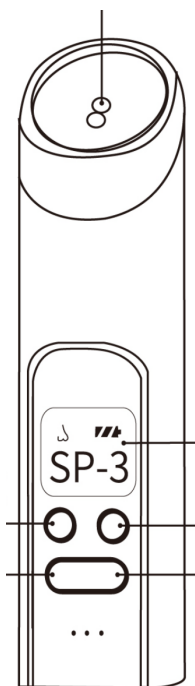
User Manual

1. Heartfelt Thank You

Thank you so much for purchasing from us and placing your trust in our brand.

GROWNSY – Parent's best helper.

GROWNSY is committed to creating modern and thoughtful baby products that let parents enjoy the parenthood. We are always passionate and devoted to making the ideal products that show love and care. Customer satisfaction is our first priority.

For any product issue, please don't hesitate to contact us.

service@grownsyofficial.com Or SMS GROWNSY to (833)524-0498
--

We will go above and beyond to resolve any customer issue within 24 hours.

Thank you again for your custom, and we truly hope you love your new nasal aspirator as much as we do.

Best Regards,


Randell

Founder of GROWNSY

2. Important Warning

- The nasal aspirator can **NOT** be turned on while charging.
- **Never allow snot collected to exceed 1/2 of the entire storage capacity! Otherwise snot will overflow.**
- **Turn off the device RIGHT AFTER snot sucking is finished, to avoid mucus back flow.**
- **If liquid enters the machine, remove the storage chamber, then run the machine for 5–10 minutes to drain the liquid that entered the motor. Wipe away stains afterwards.**
- Keep the storage chamber and the device **CLEAN** and **DRY** after use.
- To keep the battery healthy, charge the battery to 100% every two months if you do not need to use it for a long time.
- Please **ONLY** use a 5V adapter; a voltage over 5V can **NOT** charge the device.

3. Product Instruction

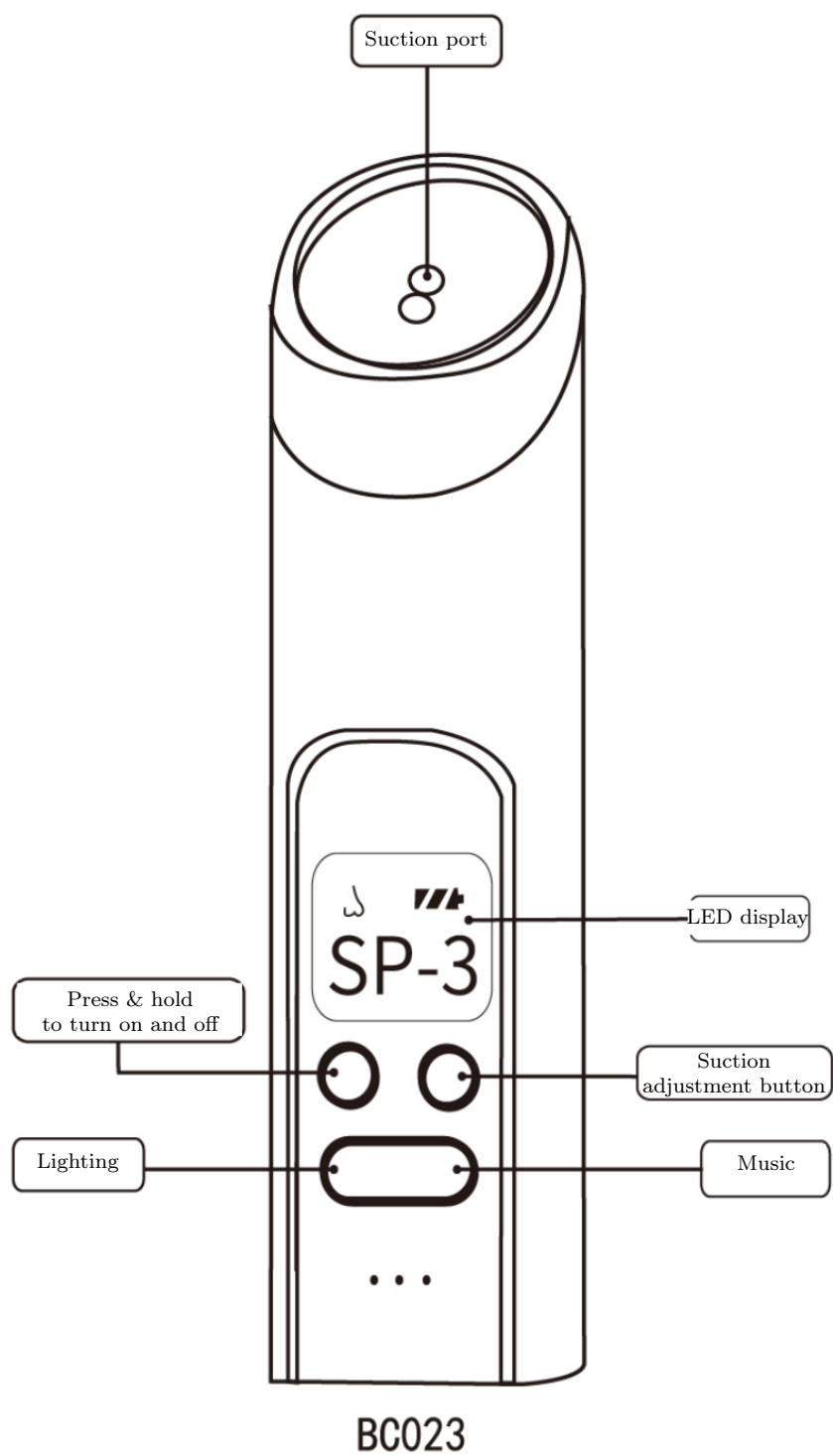4. Component Description


Funnel nozzle
For children over 2 years
old

Flat nozzle
For newborns, sucking
deep dirt


Gourd nozzle
For newborns, sucking
shallow dirt

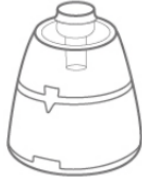Suction chamber cap

Nose suction bracket

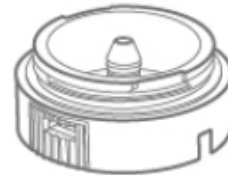Suction chamber base


Storage case


Sealant rubber ring

5. Operation of the Suction Chamber


1. Align the suction bracket with the center of the base and insert it tightly.

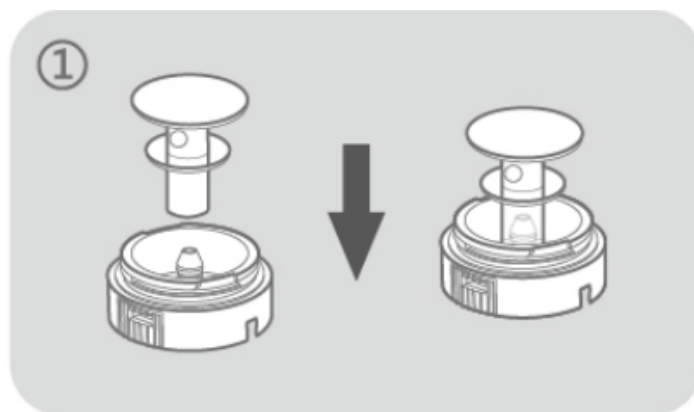
2. Fit the suction chamber onto the base and rotate it clockwise.

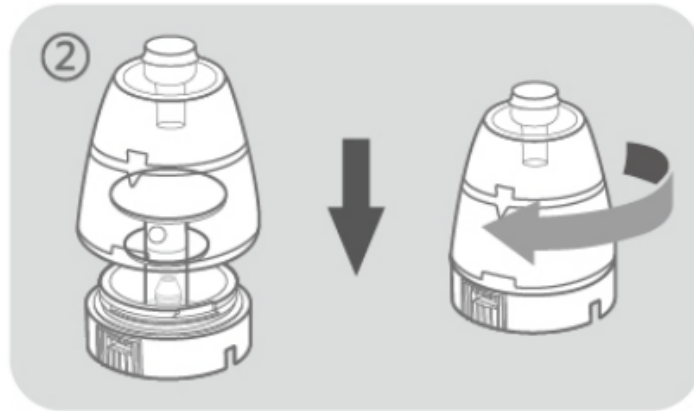3. Choose the appropriate nozzle and fit it on top of the suction chamber cap.

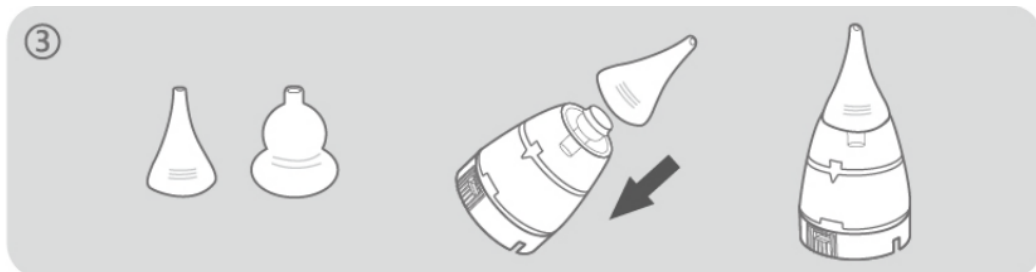
6. Connect the Suction Chamber to the Main Body

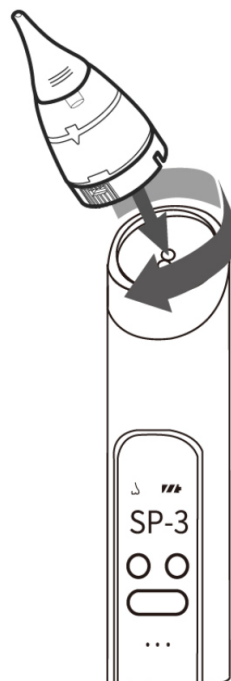Insert the assembled suction chamber set into the suction port of the main body and rotate it clockwise to ensure it fits tightly. Make sure it will NOT fall off at any angle.

Note: Suction will be weak if the sealant rubber ring falls off. Please check whether the sealant rubber ring around the suction chamber base works well. Fit it tightly before screwing the suction set into the suction port.

7. Explanation of Settings

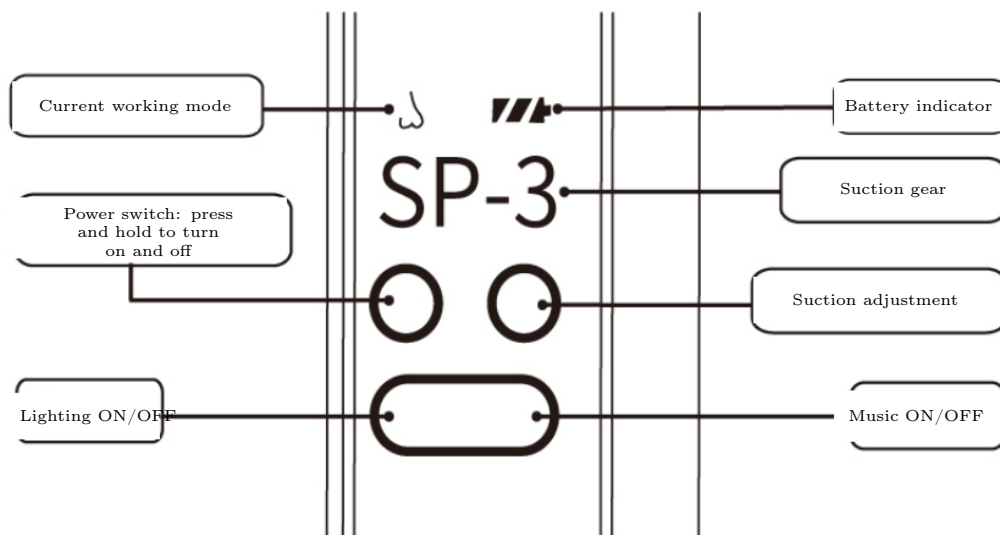
- **Power:** press and hold to turn on and off.
- **Current working mode:** nose sucking mode.
- **Suction adjustment:** there are three levels; the suction power increases from SP-1 to SP-3 level in turn. Short press to switch level.
- **Light:** the light is off by default; press the button to turn on the light, then press to turn it off.
- **Music:** the music is off by default; press once to turn on the music, press again to amplify the volume, press a third time to turn off the music.
- **Battery indicator:** will flash when there is a lack of power.

8. Note


- The nasal aspirator can not be turned on while charging; please use it while it is not being charged.
- Do not rinse or soak the main body.
- Do not charge with a voltage over 5V; such an adapter is not accepted.

Turn off the device RIGHT AFTER snot sucking is finished, to avoid mucus back flow.

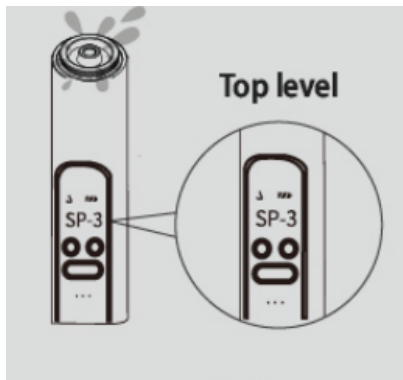If snot overflows and blocks the vent

Turn the aspirator on to suction level 3 and run it for 5–10 minutes to drain the liquid. Use a cotton swab to clean the air inlet/outlet of the main unit.

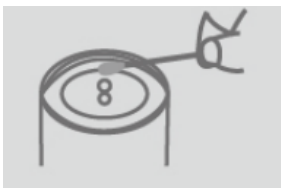Tip: Never allow the snot collected to exceed 1/2 of the entire storage capacity! Otherwise snot will overflow.

9. Precautions (Please Read Carefully Before Use)

The gourd-type silicone nozzle is suitable for sucking shallow dirt, and the flat and funnel nozzles are suitable for sucking deep dirt.

It is easier to suck out the nasal mucus if the suction tip is close to the nasal cavity, because low pressure will be formed if air cannot get in.

Please always pay attention to the nasal mucus storage chamber during use. If there is too much nasal mucus, it may overflow into the suction port, which may cause a malfunction of the product. In this condition, please turn off the machine and clean it in time. You can use a cotton swab to dry the nasal mucus from the suction port. If the suction port has sucked a lot of nasal mucus, please be sure to shake it out or turn on the device in strong suction mode to let the snot inside be sprayed out (with the suction chamber removed).

Please remove the suction chamber before cleaning. It is forbidden to immerse the product. The product is not waterproof, but the silicone nozzle can be sterilized in boiling water.

10. Specifications

Item name	Electric Baby Nasal Aspirator
Model	BC023
Power supply	3.7V
Input voltage	5V 1A
Usage time	Full charge lasts 90 min
Material	ABS + food grade PP material
Certification	FCC, CE, RoHS, REACH, WEEE, LFGB
Origin	Made in China

11. Warranty and Customer Service

GROWNSY strives for user friendly products for newborn and household items. We aim not only to create a product, but also an item that is filled with care and love from us to our customers. If you have any issues at all, please contact us directly at **service@grownsyofficial.com** and we will help you out. Or SMS GROWNSY to **(833)524-0498**.

11.1. Register Your Free Warranty

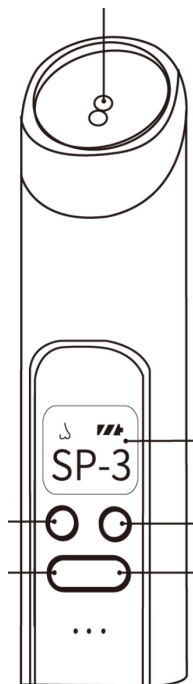
At GROWNSY we only use the highest quality product standards – but please don't forget to register for your FREE one year warranty just in case...

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