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
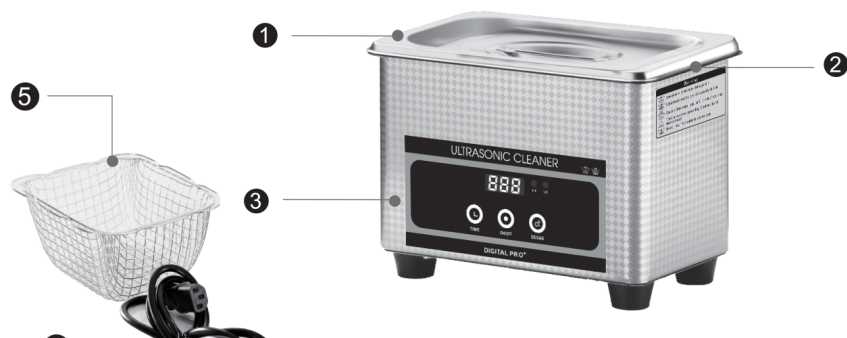
STAINLESS STEEL ULTRASONIC CLEANER

professional range · JX-008 / S / T / ST / PLUS / mechanical

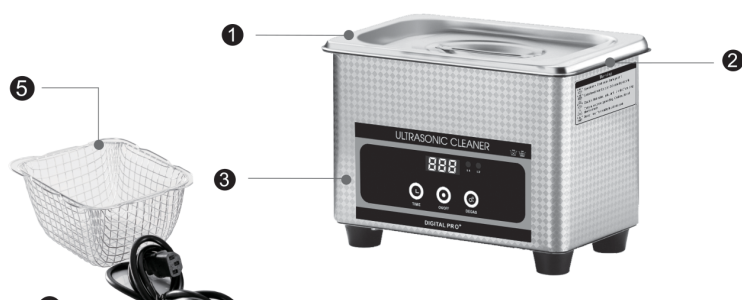User Manual

Please read this manual carefully before using the product and keep it for future reference. The pictures are for reference only; the actual product prevails.

Contents

1	Preface	3
2	The Basic Principle of Ultrasonic Cleaning	3
3	Safety Notes Before Use	3
4	Suitable For	4
5	Machine Parts	5
6	Control Panels	5
7	Operation	6
7.1	Preparation	6
7.2	Starting the ultrasound	6
7.3	Heating	7
7.4	Additional functions	7
8	Specifications	7
9	Maintenance and Troubleshooting	8
10	Disposal and Compliance	9

1. Preface

Thank you for purchasing this ultrasonic cleaner. Please read this manual carefully before use to avoid damage to the machine or danger to persons. Before plugging in the machine, make sure the power supply is within the rated range. Refitting is strictly prohibited! The control panel is eroded by organic solvents, strong acids and strong alkalis.

2. The Basic Principle of Ultrasonic Cleaning

The cleaner uses ultrasonic energy (40 kHz) in the form of sound waves to create millions of tiny microscopic bubbles in the solution, which even work their way into holes and hidden cavities, loosening dirt on all surfaces the solution touches. This action, called cavitation, occurs thousands of times per second to gently yet thoroughly scrub contamination off the item being cleaned. When you lift the item out of the cleaner, it is microscopically clean.

3. Safety Notes Before Use

1. Before connecting the power cable, make sure the power supply is within the rated range. Refitting is strictly prohibited! The control panel is eroded by organic solvents, strong acids and strong alkalis.
2. Ensure the earth wiring is well connected before starting.
3. Before starting, the power switch must be in the "0" position.
4. Do not operate with an empty tank – the ultrasonic generator will be damaged. When heating, the water level must not be below 2/3.
5. Close the lid to reduce noise, and when opening it watch out for water and steam – risk of scalding.
6. Do not relocate the machine with fluid in the tank – risk of overflow.
7. Choose the proper size and type of basket for the item to be cleaned.
8. For ultrasonic cleaning we recommend a water-soluble cleaning agent; add an appropriate amount to the tank according to the cleaning needs. Do not use flammable liquids (alcohol, petrol, etc.) – heating and long-term vibration make the liquid heat up, so beware of fire!
9. Do not use the machine in a harsh environment: where the temperature changes sharply; where humidity is too high and dew forms easily; where there is strong vibration or impact; where there is corrosive gas or dust; where water, oil or chemicals splash; where there is an explosive or flammable atmosphere.

Warning

Never start with an empty tank · the level must be above 2/3 · caution: hot liquid · earthing required · do not use flammable solutions · protect the environment when disposing.

4. Suitable For


The full range of standard ultrasonic cleaners is made of high-quality, corrosion-resistant stainless-steel plate with a long service life. It uses high-quality ultrasonic transducers and advanced bonding technology, with high electro-acoustic conversion efficiency and strong ultrasonic power output. It is equipped with an automatic constant-temperature heater in the range from room temperature to 80 °C.

Applications: electronics factories, utensils, hardware, machinery, communication, biology, medical instruments, glasses, optics, jewellery, clocks and watches, chemicals, textiles, auto parts, ships, aviation, mining and commercial cleaning.

Industry	Items and materials to be cleaned	Dirt to be removed
Semiconductors	integrated circuits, power transistors, silicon wafers, diodes, lead frames, capillaries, trays, etc.	quenching and stamping oils, polishing wax, dust particles, etc.
Electrical & electronic	picture tubes, printed circuit boards, quartz parts, electronic components, speakers, meters, LCD glass, computer parts, heads, photomasks, etc.	fingerprints, dust, cutting and stamping oils, metal filings, polishing compound, resin, etc.
Precision machinery	bearings, sewing-machine parts, typewriters, textile machines, optomechanics, gas valves, watches, cameras, metal filter elements, etc.	cutting oils, metal filings, polishing compound, fingerprints, grease, etc.
Optics	glasses, prisms, lenses, filters, glass devices, film, optical fibres, etc.	plastic, resin, paraffin, fingerprints, etc.
Hardware & machine parts	bearings, gears, shafts, tools, valves, cylinders, burners, compressors, hydraulic presses, etc.	cutting oils, metal filings, grease, polishing compound, fingerprints, etc.

Industry	Items and materials to be cleaned	Dirt to be removed
Medical instruments	medical instruments, dentures, etc.	metal filings, polishing compound, oils, stamping oils, dirt, etc.
Auto parts	piston rings, carburettors, housings, compressor bodies, electrical components, etc.	metal filings, polishing compound, oils, stamping oils, dirt, etc.
Electroplating & chemical fibre	galvanised parts, moulds, stampings; chemical-fibre nozzles and filters, etc.	chemical colloid, glue, solid residues, rust, oxide layers, dust, etc.

5. Machine Parts

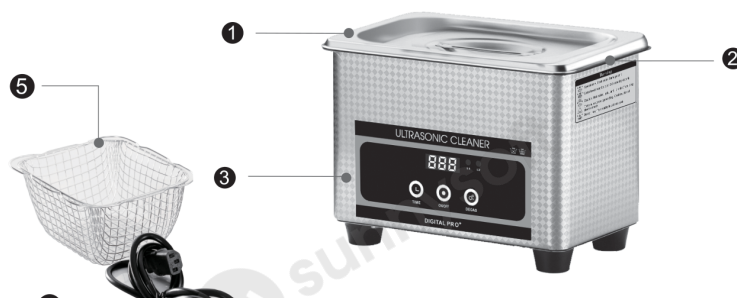- 1 Lid
- 2 Tank
- 3 Stainless-steel body (SUS)
- 4 Plug
- 5 Basket

Accessories such as cleaning baskets are shipped randomly and change according to the production batch.

Rear and housing: at the back are the drain valve, ventilation holes, carrying handle and control board.

6. Control Panels

Depending on the model, the machine has a digital panel (Digital Pro) or mechanical knobs. Functions and their availability by range:

- **ON/OFF** (ultrasonic) – turns the ultrasound on and off; the ultrasonic indicator lights up.
- **TIME** or **TIME +/TIME –** – sets the cleaning time (fixed settings or continuous, depending on the range); the display shows a countdown.

- **HEAT / HEATER ON/OFF** and **TEMP +/TEMP –** – turns the heater on and sets the temperature (from room temperature to 80 °C); the heater indicator lights up.
- **DEGAS** – degasses the liquid.
- **SWEEP** – frequency sweeping (ST range) for more even cleaning.
- **SEMIWAVE** – half-wave (half the ultrasonic power) for cleaning delicate objects.
- **PLUS** – power boost (PLUS range) for stubborn stains.
- Mechanical range: time knob (0–30 min), power knob and heater switch.

7. Operation

Preparation

1. Add an appropriate amount of cleaning liquid to the tank. With heating on, the level must not fall below 2/3 of the tank height. Do not use flammable liquids (alcohol, petrol, etc.) – heating and vibration make them heat up, so beware of fire! Strong acids and alkalis are strictly forbidden – risk of chemical reaction and corrosion of the tank!
2. Plug the machine into an earthed socket and turn on the switch at the back ("0" = off, "I" = on).

Starting the ultrasound

JX-008 / JX-008T Use **TIME** to select the cleaning time. There are 18 settings: 90 / 120 / 150 / 180 / 210 / 240 / 270 / 300 / 330 / 360 / 540 / 720 / 900 s and 18 / 21 / 24 / 27 / 30 min. Press **ON/OFF** to start; press again to pause.

JX-009 / JX-010T Set the cleaning time; there are 7 values to choose from: 60 / 180 / 300 / 600 / 900 s, 30 min and 60 min. After setting, press **ON/OFF** – the ultrasonic indicator lights up and the machine works until the time is up. Press **ON/OFF** again to pause.

JX-S / JX-PLUS / JX-T / JX-ST After switching on, a random time is shown. Use **TIME +** to increase the set time by 1 minute and **TIME –** to decrease it by 1 minute. After setting, press **ON/OFF** – the ultrasonic indicator turns red and the machine works according to the set parameters; the timer counts down and stops automatically at "00:00".

Mechanical range Turn the knob (ULTRASONIC) clockwise to select a time of 0–30 min. When the indicator lights up and a "ZIZI" sound is heard, the ultrasound is working.

Heating

JX-009 / JX-010T The heater indicator lights up after pressing **HEAT**; the heater switches off automatically when 65 °C is reached. Press again to pause.

JX-S / JX-PLUS / JX-T / JX-ST After switching on, a random temperature is shown, which can be set. **TEMP +** increases the temperature by 1 °C, **TEMP -** decreases it by 1 °C. After setting, press **ON/OFF** – the machine starts heating and, when the set temperature is reached, the heater indicator goes out.

Mechanical range If heating is needed, start it (HEATING) and set the required temperature; 40–60 °C is usually best. Heating is optional according to need.

Additional functions

DEGAS on pressing, the machine works for 9 seconds and stops for 6 seconds (degassing). Range JX-T / 008T / ST / 009, 010T.

SWEEP on pressing, the machine sweeps the frequency within the action range; the liquid forms a fine backflow and the dirt is removed from the workpiece faster. Range JX-ST.

SEMIWAVE turns on the half-wave – the ultrasonic power is halved; for delicate objects. Range JX-T.

PLUS on pressing during cleaning, increases the power (without affecting the machine) for stubborn stains. Range JX-PLUS.

Note

When finished, press the ultrasonic **ON/OFF** to stop the ultrasound (the indicator goes out) and **HEATER ON/OFF** to stop the heater. Then unplug the machine, pour out the liquid and wipe the tank and the machine with a cloth.

8. Specifications

The dimensions are measured by hand; the actual product prevails. The dimensions given do not include the drain valve.

Small models (JX-008 range):

Model	Vol. (l)	Tank (mm)	Unit (mm)	Ultras. (W)	Heat (W)	Time	Temp. (°C)
JX-008	0.8	150×85×65	178×110×130	35	–	18 settings	auto 65
JX-008T	0.8	150×85×65	178×110×130	35	–	18 settings	auto 65
JX-009	1.3	150×135×65	175×160×175	60	40	18 settings	RT-80
JX-009S	1.3	150×135×65	175×160×175	60	40	7 settings	RT-80
JX-010S	2	150×135×100	175×160×210	100	–	–	RT-80
JX-010T	2	150×135×100	175×160×210	100	60	7 settings	auto 65

Large models (JX-020 to JX-100 – mechanical, S, T, ST):

Model	Vol. (l)	Tank (mm)	Unit (mm)	Heat (W)	Ultras. (W)	Temp. (°C)
JX-010	2	150×135×100	175×160×210	100	60	RT-80
JX-020	3.2	240×135×100	265×165×220	100	120	RT-80
JX-030	4.5	300×150×100	325×180×225	200	180	RT-80
JX-031	6.5	300×150×150	325×180×280	200	180	RT-80
JX-040	10	300×240×150	325×265×280	200	240	RT-80
JX-060	15	330×300×150	360×325×285	300	360	RT-80
JX-080	22	500×300×150	530×325×285	500	480	RT-80
JX-100	30	500×300×200	530×325×325	500	600	RT-80

Frequency of all models 40 kHz. Key: "settings" = number of cleaning-time settings; "RT-80" = from room temperature to 80 °C; "auto 65" = automatic heating to 65 °C. The T and ST ranges have selectable (half / full) ultrasonic power; the table shows the full power.

The PLUS range additionally has a boosted PLUS power (e.g. up to 180 W on the JX-020PLUS).

9. Maintenance and Troubleshooting

Regularly blow the dust off the electrical parts with compressed air and clean the tank at least once a week.

Troubleshooting:

Problem	Possible cause	Solution
No ultrasound	power supply not connected; fuse blown; short circuit in cable or transducer; PCB board faulty	check and connect the power supply; replace the fuse with the same type; replace the cable; replace the faulty part; otherwise contact service
No heating	heater switch faulty; fuse blown; heating pad faulty; temperature set incorrectly	check the heater wiring; replace the fuse and the heating pad; the set temperature must be 3 °C higher than the actual one; otherwise contact service
Poor cleaning	level too high / too low; unsuitable temperature; unsuitable liquid	adjust the level and temperature; after cooling, replace with a suitable liquid; otherwise contact service
Timer failure	timer faulty	replace the timer or the digital panel; otherwise contact service

Problem	Possible cause	Solution
Temperature-control failure	thermostat loose or broken; display out of control	tighten or replace the thermostat; replace the faulty part; otherwise contact service
Current leakage	incorrect earthing; live or neutral wire loose	check the earthing and internal wiring; contact service

10. Disposal and Compliance

 The crossed-out wheeled-bin symbol means that the product does not belong in normal household waste (WEEE under Directive 2012/19/EU). At the end of its service life, hand it in at a collection point for electronic waste. In this way you protect human health and the environment.

 This product complies with the applicable requirements of the European Union (CE) and the RoHS Directive on the restriction of hazardous substances.

Manufacturer: Shaoguan Skymen Ultrasonic Technology Limited

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